

EYA2 mAb

Catalog Number: 26096

Gene Symbol: EYA2

Description: Anti- EYA2 Mouse Monoclonal Antibody

Background: EYA2 may be involved in development of the eye. It promotes efficient DNA repair by dephosphorylating H2AX, promoting the recruitment of DNA repair complexes containing MDC1. Its function as histone phosphatase probably explains its role in transcription regulation during organogenesis. EYA2 coactivates SIX1 and seems to coactivate SIX2, SIX4 and SIX5. It seems to be involved in myogenesis together with SIX1 and DACH2.

Immunogen: Human EYA2 protein

Tested applications: ELISA, WB, IHC

Recommended dilutions:

ELISA: 1:1000-1:5000

WB: 1:100-1:1000

IHC: 1:50-1:100

Concentration: 0.8 mg/ml

Host: Mouse

Clonality: Monoclonal

Purity: Purified from ascites

Format: Liquid

Storage buffer:

Preservative: no


Constituents: PBS (without Mg²⁺ and Ca²⁺), pH 7.4, 150 mM NaCl, 50% glycerol

Species Reactivity: Anti-EYA2 antibody recognizes EYA2 of vertebrates.

Storage Conditions: Store at -20 °C. Avoid freeze / thaw cycles

Western blot:

EYA2 protein

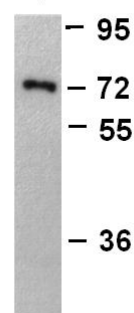


WB: Anti-EYA2 mAb

Western blot analysis of recombinant full length EYA2 protein. Purified His-tagged EYA2 protein (Cat. #10195) was blotted with anti-EYA2 monoclonal antibody (Cat # 26096).

For research use only. Not for diagnostic or therapeutic applications.

Immunohistochemistry:

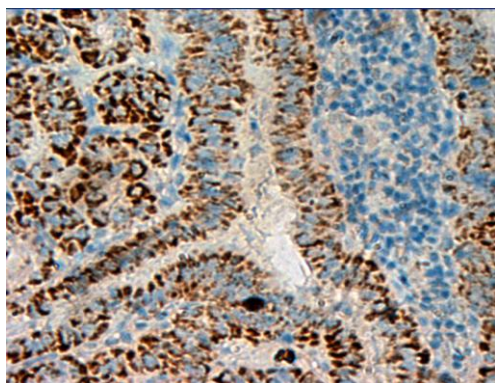

Immunohistochemical analysis of paraffin-embedded thyroid cancer tissue with anti-EYA2 monoclonal antibody (Cat. #26096).

Tissue samples were fixed with formaldehyde and blocked with 1% serum for 15 min at 37 °C. Antigen retrieval was by heat mediation in citrate buffer (pH6). Samples were then incubated with primary antibody (1:50) overnight at 4 °C. A HRP-conjugated Goat anti-mouse IgG (dilution 1:50) was used as secondary antibody.

For research use only. Not for diagnostic or therapeutic applications.