

Product Name: Sar1a Protein H79G mutant

Synonyms: SAR1 homolog A, SAR1, Sara, SARA1, Masra2

Catalogue Number: 10115

Source: Human, recombinant full length, His₆-tag

Expression Host: E. coli

Molecular Weight: 22 kDa

Purity: >95% by SDS-PAGE

Introduction: Small GTPases are a super-family of cellular signaling regulators. Sar1 is a member of these small GTPases. Sar1 functions as a molecular switch to control protein–protein and protein–lipid interactions that direct vesicle budding from the ER.

Amino Acid Sequence (1-198, H79G)

MSFIFEWIYNGFSSVLQFLGLYKKSGLVFLGLDNAGK
TLLHMLKDDRLGQHVPTLHPTSEELTIAGMTFTTDFL
GGGEQARRVWKNYLPAINGI VFLVDCADHSRLVESKVE
LNALMTDETISNVPILILGNKIDRTDAISEEKLREIFG
LYGQTTGKGNVTLKELNARPMEVFMCSVLKRQGYGEGF
RWLSQYID

Properties

Physical Appearance (form): Lyophilized Powder from 20 mM Tris-HCl, pH 7.4, containing 100 mM NaCl.

Physical Appearance (color): White


Solubility (Color): colorless

Solubility (Turbidity): clear, 10 mg/mL H₂O

Storage: -80°C

Preparation Instructions

Centrifuge the vial before open the cap and reconstitute in water. Adding of 10 mM β-mercaptoethanol or 1 mM DTT into the solution to protect the protein is recommended and using of non-ionic detergents such as n-Dodecyl β-D-maltoside (DoDM) or polyethylene detergents (e.g. C₁₂E₁₀) also help to stabilize the protein. Avoid repeated freezing and thawing after reconstitution.

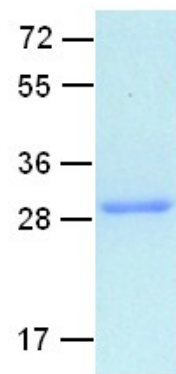

References

1. He, H. et al., Gene Expr. 10: 231-242, 2002.
2. Jones, B. et al., Nature Genet. 34: 29-31, 2003.
3. South, S. T. et al., J. Cell Biol. 149: 1345-1359, 2000.

For research use only. Not for diagnostic or therapeutic applications.